





Lyo-Ready 1-Step RT-qPCR Mixes

Lyophilization-compatible RT-qPCR master mixes that are ideal for high-throughput and point-of-care diagnostic assays.

Suitable Applications

-  Environmental
-  Human Diagnostics
-  Vet Health
-  DNA Barcoding

Lyo-Ready 1-Step RT-qPCR master mixes are glycerol-free formulations that contain Taq polymerase, reaction buffer, dNTP, MgCl₂, lyo-exipients and separate reverse transcriptase and dilution buffer. They are ideal for developing room temperature, lyophilized, multiplex RT-qPCR diagnostic tests and are suited for high-throughput, automated platforms requiring milli-second rehydration times.

Optimized 1-Step RT-qPCR Master Mixes for Lyophilization


- Provided glycerol-free and pre-formulated with a specialized blend of lyo-exipients
- Virus mix (Cat #MDX062) is specifically optimized for amplification of RNA and DNA viruses
- Suited for multiplex-assays and low-copy number targets
- Compatible with lyophilization into beads or cakes
- Ideal for developing fast, highly reproducible assays requiring ambient temperature transport and storage

PRODUCT	CAT NO.	VOLUME	REACTIONS
Lyo-Ready 1-Step RT-qPCR Mix	MDX024	10 mL	1,000 Rxn
		100 mL	10,000 Rxn
Lyo-Ready 1-Step RT-qPCR Virus Mix	MDX062	10 mL	1,000 Rxn
		100 mL	10,000 Rxn

Product Highlights


High efficiency and sensitivity from both DNA and RNA templates

Amplification of both Rotavirus A (dsRNA) and Adenovirus (dsDNA) in a single multiplexed RT-qPCR assay, from inactivated crude viral lysates. The result illustrates that both viruses were amplified with high efficiency and sensitivity, demonstrating the ability of Lyo-Ready 1-Step RT-qPCR Virus Mix (Cat #MDX062) to detect low-copy number RNA and DNA targets simultaneously from a single sample.

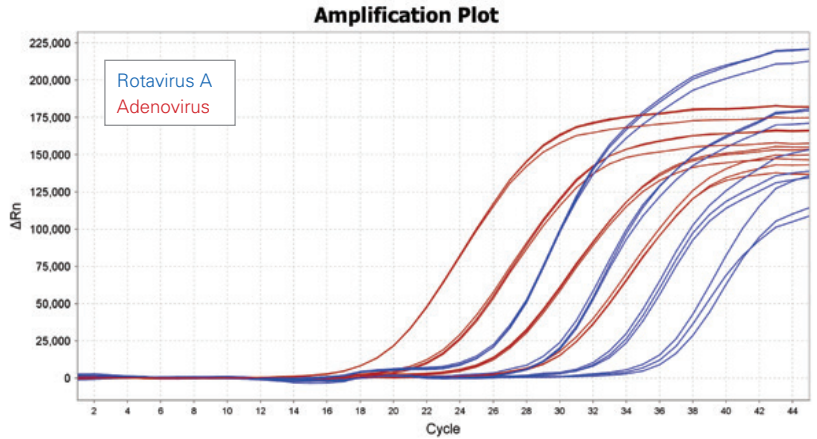




Inhibitor tolerant enabling amplification directly from blood samples

γ -actin was amplified from mouse RNA in the presence of 2% whole blood using Lyo-Ready 1-Step RT-qPCR Mix (purple traces) and a competitor mix (green traces). The results demonstrate the ability of Lyo-Ready 1-Step RT-qPCR Mix (Cat #MDX024) to amplify from difficult samples such as blood, in contrast to the competitor mix which fail, even with higher concentrations of the template RNA.

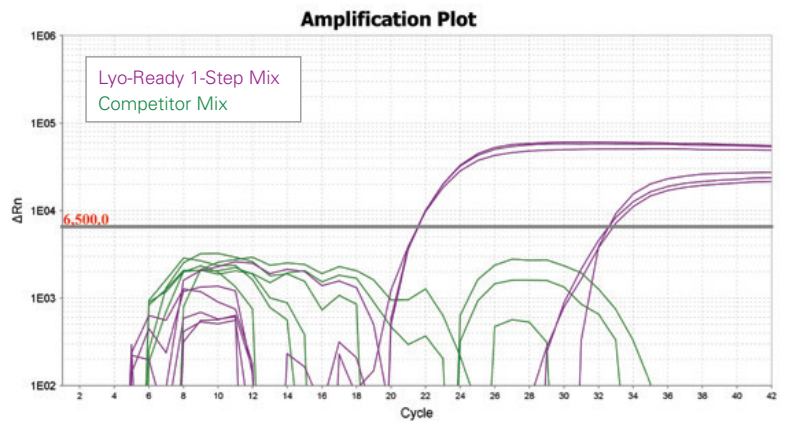

Retention of reverse transcriptase activity following lyophilization

Primers and probes were added to Lyo-Ready 1-Step RT-qPCR Mix (Cat #MDX024), dried down and immediately rehydrated (purple) and tested against a freshly prepared wet mix (orange), in a triplex probe RT-qPCR assay on mouse RNA. The results demonstrate the same reverse transcriptase enzyme activity before and following lyophilization.

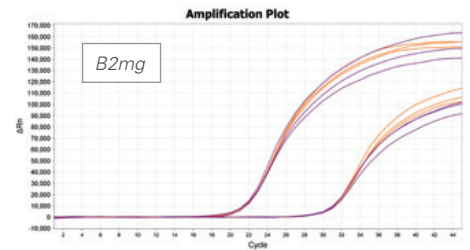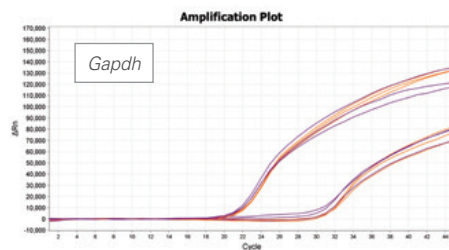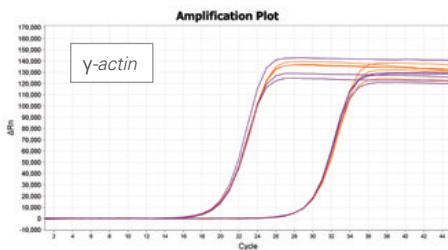

Lyophilised | Wet Mix

Ordering information:

USA
E: info@meridianlifescience.com
Toll free: +1 800 327 6299

Germany
E: info.de@bioline.com
Tel: +49 (0)3371 60222 00



美国迈迪安生命科学公司
电子邮件: vivian.li@meridianlifescience.com
电话: +86-159-1103-0750

UK
E: info.uk@bioline.com
Tel: +44 (0)20 8830 5300

France
E: info.fr@bioline.com
Tel: +33 (0)1 42 56 04 40

Australia
E: info.au@bioline.com
Tel: +61 (0)2 9209 4180

meridian BIOSCIENCE™
LIFE DISCOVERED. LIFE DIAGNOSED.

Connect with us:  

www.MeridianBioscience.com/Lifescience

ISO 13485 Certified

01/19